

Nano-Secondary[®] alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor[®] 488 [CTK0105, CTK0106]



Product code: sms2bAF488-1

Properties

Description	Monovalent, recombinant secondary single domain antibodies to mouse IgG2b: Mixture of 2 alpaca monoclonal Nanobodies, Fc-specific, Alexa Fluor [®] 488 conjugated
Product Type	Nano-Secondary [®] , secondary Nanobody (VHH)
Format	Alpaca single domain antibody, monovalent
Host	Alpaca-derived, recombinantly produced in bacteria
Target / Specificity	Fc-fragment of mouse IgG2b
Cross-reactivity	No cross-reactivity to goat, guinea pig, human, macaque, rabbit, rat, and sheep serum and to IgG1, IgG2a, IgG2c, and IgG3
Immunogen	Purified mouse IgG2b
Clonality	Biclonal: mixture of 2 monoclonal Nanobodies
Clones	CTK0105, CTK0106 (VHH0275, VHH0288)
Affinity (Kd)	CTK0105: Kd = 5.0 nM, CTK0106: Kd = 0.2 nM
Conjugate	Site-directed conjugation to Alexa Fluor 488
Excitation/ Emission	Excitation max: 490 nm, Emission max: 525 nm
Degree of Labeling (DOL)	2 fluorophores per Nanobody
Synonyms	Alpaca single domain antibody, VHH, Nanobody, binding domain of single domain antibody, Nano-antibody
Validation	Application validated for immunofluorescence and Western blotting Determination of cross-reactivity, sequence, affinity, melting point, and degree of labeling (DOL)
Purity	Recombinantly expressed and purified
Form	Buffered aqueous solution
Concentration	0.5 mg/mL
Storage buffer	10 mM HEPES pH 7.0, 500 mM NaCl, 5 mM EDTA Preservative: 0.09 % sodium azide, safety datasheet (SDS): sodium azide
Storage instructions	Shipped at ambient temperature. Store at -20°C/-4°F. Avoid freeze-thaw cycles. Aliquot upon arrival. Protect from light. Stable for 6 months.
Size	10 µL; 100 µL
RRID	AB_2827581

Nano-Secondary[®] alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor[®] 488 [CTK0105, CTK0106]




Product code: sms2bAF488-1

Applications	<p>Spin down material before use!</p> <p>IF/ICC: recommended starting dilution 1:1,000 (e.g. PBS supplemented with 4% BSA)</p> <p>Western blot: recommended starting dilution 1:1,000 (e.g. PBS supplemented with 0.075% Tween-20 and 5% skimmed milk)</p> <p><i>Note:</i> The optimal dilution depends on the application and should be determined by the user. A titration from range from 1:250 up to 1:2,000 is recommended.</p> <p><i>Note:</i> Image acquisition time may have to be optimized.</p>
---------------------	-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------


Tested Applications

Immunofluorescence

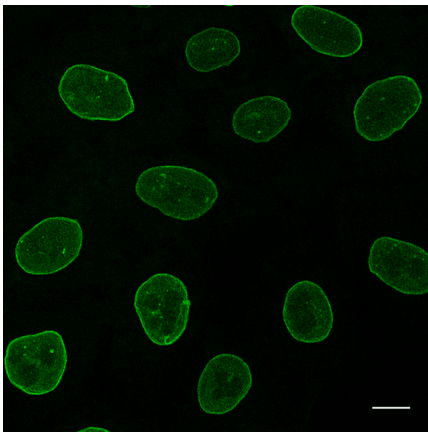

Immunostaining of Lamin in HeLa cells with mouse anti-Lamin A/C antibody and Nano-Secondary[®] alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor[®] 488 [CTK0105, CTK0106] 1:1,000 (green). Scale bar, 10 μ m.

Western Blot

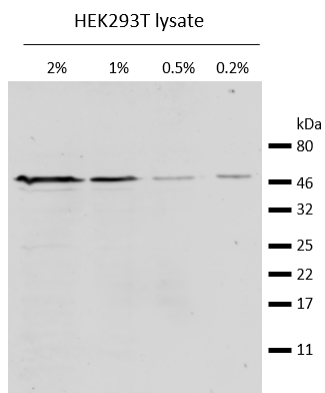

Western blot analysis of endogenous β -Tubulin in HEK293T cell lysate. Detection with mouse anti- β -Tubulin antibody and Nano-Secondary[®] alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor[®] 488 [CTK0105, CTK0106] 1:1,000.

Nano-Secondary® alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor® 488 [CTK0105, CTK0106]



Product code: sms2bAF488-1

Product overview and related products

Nano-Secondaries (Alexa Fluor® conjugated)	Product code
Nano-Secondary® alpaca anti-human IgG/anti-rabbit IgG, recombinant VHH, Alexa Fluor® 488 [CTK0101, CTK0102]	srbAF488-1-10; -100
Nano-Secondary® alpaca anti-human IgG/anti-rabbit IgG, recombinant VHH, Alexa Fluor® 568 [CTK0101, CTK0102]	srbAF568-1-10; -100
Nano-Secondary® alpaca anti-human IgG/anti-rabbit IgG, recombinant VHH, Alexa Fluor® 647 [CTK0101, CTK0102]	srbAF647-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 488 [CTK0103, CTK0104]	sms1AF488-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 568 [CTK0103, CTK0104]	sms1AF568-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, Alexa Fluor® 647 [CTK0103, CTK0104]	sms1AF647-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor® 488 [CTK0105, CTK0106]	sms2bAF488-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor® 568 [CTK0105, CTK0106]	sms2bAF568-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor® 647 [CTK0105, CTK0106]	sms2bAF647-1-10; -100
Nano-Secondary® alpaca anti-mouse IgG3, recombinant VHH, Alexa Fluor® 647 [CTK0107]	sms3AF647-1-10; -100
Nano-Secondaries (unconjugated)	
Nano-Secondary® alpaca anti-human IgG/anti-rabbit IgG, recombinant VHH, for 2x Cys conjugation [CTK0101]	srbGCys2-1-500
Nano-Secondary® alpaca anti-rabbit IgG, recombinant VHH, for 2x Cys conjugation [CTK0102]	srbGCys2-2-500
Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, for 2x Cys conjugation [CTK0103]	smsG1Cys2-1-500
Nano-Secondary® alpaca anti-mouse IgG1, recombinant VHH, for 2x Cys conjugation [CTK0104]	smsG1Cys2-2-500
Nano-Secondary® alpaca anti-mouse IgG2b, recombinant VHH, for 2x Cys conjugation [CTK0105]	smsG2bCys2-1-500
Nano-Secondary® alpaca anti-mouse IgG2b, recombinant VHH, for 2x Cys conjugation [CTK0106]	smsG2bCys2-2-500

For product details, information, and ordering visit www.chromotek.com.

Nano-Secondary® alpaca anti-mouse IgG2b, recombinant VHH, Alexa Fluor® 488 [CTK0105, CTK0106]



Product code: sms2bAF488-1

Contact

support@chromotek.com

ChromoTek GmbH
Am Klopferspitz 19
82152 Planegg-Martinsried
Germany
phone: +49 89 124 148 80
fax: +49 89 124 148 811

ChromoTek Inc.
62-64 Enter Lane
Islandia, NY 11749
USA
phone: 631 501 1058
fax: 631 501 1060

Disclaimer

Only for research applications, not for diagnostic or therapeutic use!

ChromoTek and GFP-Trap, RFP-Trap, Myc-Trap, Spot-Trap, Spot-Tag, Spot-Label, Spot-Cap, Nano-Secondary, F2H Kit, and Chromobody are registered trademarks of ChromoTek GmbH, part of Proteintech Group. Nano-CaptureLigand and V5-Trap are trademarks of ChromoTek GmbH, part of Proteintech Group. Nanobody is a registered trademark of Ablynx, a Sanofi company. Alexa Fluor is a registered trademark of Life Technologies Corporation, a part of Thermo Fisher Scientific Inc. Dynabeads is a trademark of Life Technologies AS, a part of Thermo Fisher Scientific Inc. SNAP-tag is a registered trademark and CLIP-tag is a trademark of New England Biolabs, Inc. Octet is a registered trademark of FortéBio, a Sartorius brand. Other suppliers' products may be trademarks or registered trademarks of the corresponding supplier each. Statements on other suppliers' products are given according to our best knowledge.